


Nivvi long-term support guide

A practical guide for compatible Bluetooth heart-rate monitoring, troubleshooting and release validation.

Private review draft. Replace the operator and support details before public release. Nivvi is not a medical device and does not replace a clinical plan or emergency care.

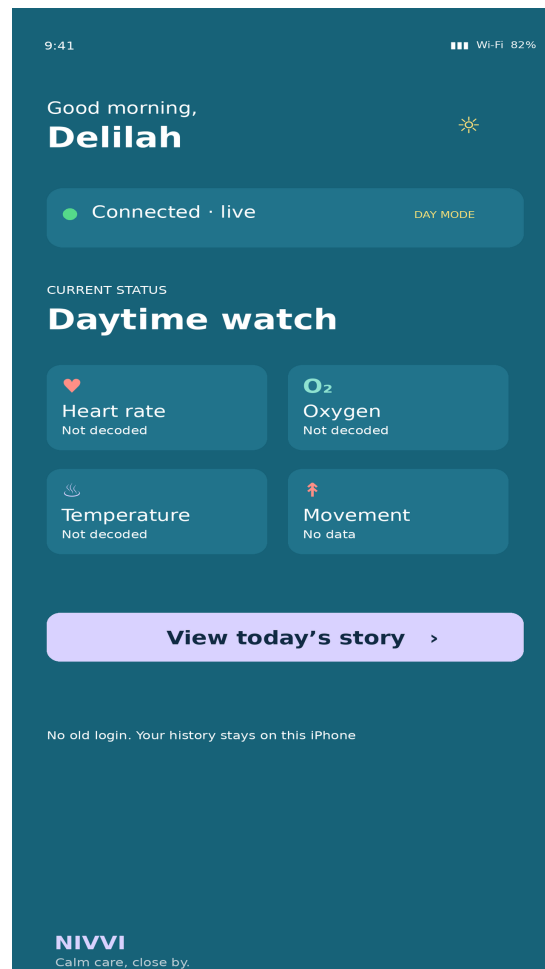

Compatibility at a glance

Standard Bluetooth Heart Rate Service devices using service 180D and characteristic 2A37 are the primary supported format. Battery percentage is shown when Battery Service 180F / 2A19 is available. Other adapter formats remain experimental and require independent verification.

Daily use

Connect and stay connected

Favourite a device in Device so it stays at the top of the scan list. Leave the session connected until you deliberately tap Disconnect. Nivvi attempts reconnection after signal loss, but range, battery, force-quitting, Bluetooth permissions and iOS restrictions can interrupt it.

History and events

History retains 30 calendar days on the iPhone. Readings are sampled every 30 seconds for charts and export, while every valid incoming standard reading is used for alarm evaluation. Events include connection changes, alarms, alarm-clear events and parent notes.

No fresh measurements

Check that the device is worn correctly, has skin contact, is charged and is not connected to another Bluetooth app. A connected status alone does not prove that a reading is available or accurate.

Alerts and sharing

Care-plan limits

Enter the low and high limits supplied by the child's care team. Low means strictly below the low limit; high means strictly above the high limit. The reading must remain outside the selected range for the configured duration. A fresh in-range reading clears the active alarm. Test the siren and lock-screen notification before relying on them.

iOS delivery limits

Silent mode, Focus, notification permissions, volume, background delivery and Critical Alerts entitlement decisions are controlled by iOS. Nivvi cannot guarantee audible delivery in every state.

Family sharing

Readings, events and notes stay on the iPhone by default. Live remote sharing should only be enabled after secure accounts, consent, revocation, stale-data indicators, access logging and monitored support are implemented.

Release validation checklist

Before a paid monitoring release, physically test each item below on every device model named in the store listing:

Area	Pass condition
Discovery	Device appears through the supported Bluetooth service path.
Reading	Valid values appear and invalid/contact-loss packets are rejected.
Alarm	Low and high alarms trigger only after the configured dwell and clear on a fresh in-range reading.
Locked phone	Notification and permitted sound behaviour are recorded under lock.
Reconnection	Signal loss is logged and the saved favourite reconnects.
History	24-hour day selection, timestamps, event filter, chart zoom/selection, export and deletion all work.
Safety	Store copy, privacy notice, support instructions and review notes match tested behaviour.

Never change a child's condition to test an alarm. Use a controlled test peripheral or recorded fixture and document the device, firmware, iOS version and result.